



vison Documentation

Release 0.2+28.gd877ee4.dirty

Ruyman Azzollini

Sep 21, 2017

CONTENTS

1	README	3
2	Installation	5
2.1	Installation	5
2.2	Dependencies	7
3	Pipeline Core	9
3.1	Pipeline	9
4	Data Model	11
4.1	Data Model	11
5	Analysis (Shared)	15
5.1	Analysis (Shared)	15
6	Charge Injection Tools	19
6.1	Charge Injection Tools	19
7	“Flat” Acq. Analysis Tools	21
7.1	“Flat” Acq. Analysis Tools	21
8	Image	23
8.1	Image Analysis	23
9	Monitoring (“Eyegore”)	25
9.1	Monitoring (“Eyegore”)	25
10	OGSE	27
10.1	OGSE Tools	27
11	Point-Source Analysis	29
11.1	Point-Source Analysis	29
12	Test Scripts	35
12.1	Charge Injection Scripts	35
12.2	Dark Scripts	36
12.3	Flat-Illumination Scripts	37
12.4	Other Test Scripts	40
12.5	Point-Source Scripts	40
12.6	Trap-Pumping Scripts	44

13 Scripts	45
13.1 Scripts	45
14 Support Code	47
14.1 Support Code	47
15 Indices and tables	51
Python Module Index	53

Contents:

README

vison Euclid VIS Ground Calibration Pipeline

Author Ruyman Azzollini

Contact r.azzollini_at_ucl.ac.uk

issue 0.2+28.gd877ee4.dirty

version 0.2+28.gd877ee4.dirty

date Sep 21, 2017

This Python package “vison” is the pipeline that will be used at MSSL for ground calibration of the VIS detection chains, including one ROE, one RPSU and three CCDs.

INSTALLATION

The package is distributed via github. The repository is hosted at:

<https://github.com/ruymanengithub/vison>

Detailed instructions:

2.1 Installation

2.1.1 Cloning *vison* from the repository using *git*

If you don't have *git* installed in your system, please follow this [link](#) first.

Here we will follow these [instructions](#) to clone the repository to your own computer. Follow the link for instructions in other operative systems.

Step-by-step:

- Go to <https://github.com/ruymanengithub/vison>.
- Click on the green "Clone or download" button.
- In the Clone with HTTPs section, click to copy the clone URL for the repository.
- Open a Terminal.
- Change the current working directory to the location where you want the cloned directory to be made.
- Type `git clone`, and then paste the URL you copied in Step 1.


```
~$ git clone https://github.com/ruymanengithub/vison
```

- Press Enter. Your local clone will be created.

2.1.2 Installation

We recommend installing the code through a *conda* environment, with a specific list of packages, so you can be sure you have all the needed dependencies.

First, if you don't have *conda* already installed in your system already, follow the instructions in this [link](#).

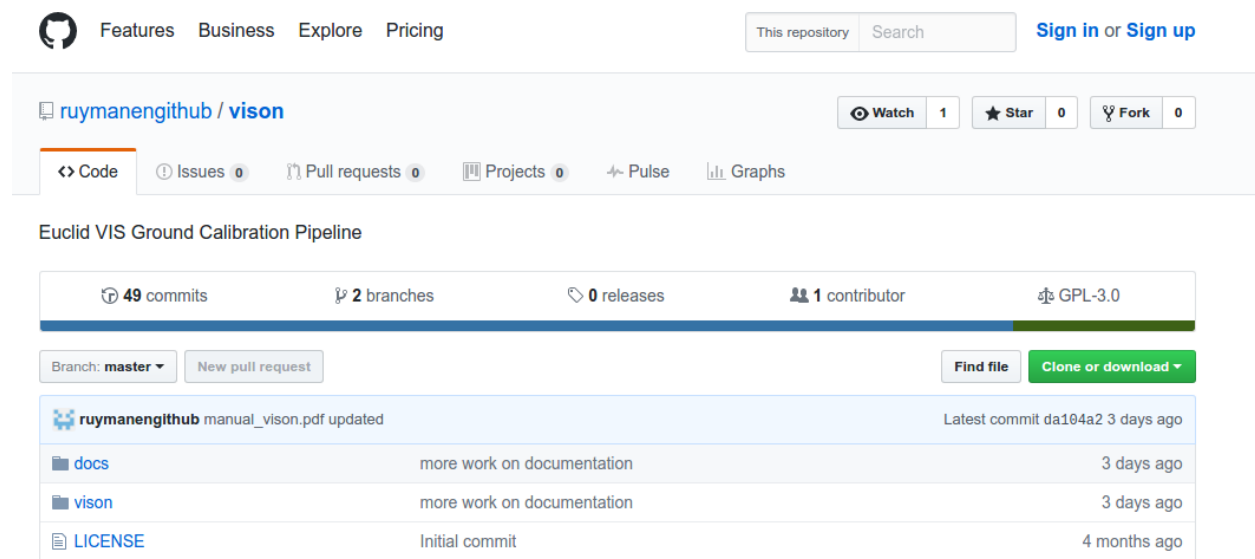

GitHub repository page for **ruymanengithub / vison**. The page displays repository statistics: 49 commits, 2 branches, 0 releases, 1 contributor, and GPL-3.0 license. It also shows a list of files: **docs** (more work on documentation, 3 days ago), **vison** (more work on documentation, 3 days ago), and **LICENSE** (Initial commit, 4 months ago).

Installing conda and creating *vison* environment

Once you have successfully installed conda, we will create an environment that will allow you to install the pipeline and meet all its dependencies.

Step-by-Step:

- change directory to your copy of the vison repository:

```
~$ cd vison
```

- Under the 'conda' sub-folder, you will find two text files:

```
~$ cd conda
~$ ls
env-conda_vison.txt  env-pip_vison.txt
```

- Then execute the following command to create a new conda environment, *vison*:

```
~$ conda create -n vison --file env-conda_vison.txt
```

- When prompted, type "y" and return to install the listed packages.
- Activate the new environment

```
~$ source activate vison
```

- Install the packages that are accessed via *pip*, within the conda environment:

```
~$ pip install -r env-pip_vison.txt
```

Installing *vison*

Finally, to install the *vison* pipeline itself, we will go back to the folder we downloaded from the github repository:

```
~$ cd ../
~$ ls
conda  docs  LICENSE  manual_vison.pdf  README.md  setup.cfg  setup_distutils.py
↪setup.py  vison
```

Then do the actual installation, via:

```
~$ python setup.py install
```

Now the vison package will be accessible from anywhere in your system, whenever you start python from within the *vison* conda environment. For example:

- open a new terminal and go to your home directory

```
~$ cd
```

- activate the vison environment:

```
~$ source activate vison
```

- start the python interpreter and import vison:

```
~$ source activate vison
~$ python
>>> import vison
>>> dir(vison)
['Eyegore', 'FlatFielding', 'Pipe', 'Report', '__all__', '__builtins__', '__doc__'
↪, '__file__',
'__name__', '__package__', '__path__', '__version__', 'analysis', 'data',
↪ 'datamodel',
'eyegore', 'pipe', 'point', 'support']
```

2.2 Dependencies

Instructions to acquire a copy of the “conda” environment that provides all dependencies is included in the package. See [Installation](#) instructions for details.

PIPELINE CORE

3.1 Pipeline

3.1.1 master.py

This is the main script that will orchestrate the analysis of Euclid-VIS FM Ground Calibration Campaign.

The functions of this module are:

- Take inputs as to what data is to be analyzed, and what analysis scripts are to be run on it.
- Set the variables necessary to process this batch of FM calib. data.
- Start a log of actions to keep track of what is being done.
- Provide inputs to scripts, execute the analysis scripts and report location of analysis results.

Some Guidelines for Development:

- Input data is “sacred”: read-only.
- Each execution of Master must have associated a unique ANALYSIS-ID.
- All the Analysis must be divided in TASKS. TASKS can have SUB-TASKS.
- All data for each TASK must be under a single directory (TBC).
- All results from the execution of FMmaster must be under a single directory with subdirectories for each TASK run.
- **A subfolder of this root directory will contain the logging information:** inputs, outputs, analysis results locations.

Created on Wed Jul 27 12:16:40 2016

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

class vison.pipe.master.**Pipe** (*inputdict, dolog=True*)

Master Class of FM-analysis

BIAS01 = <module ‘vison.dark.BIAS01’ from ‘/home/raf/SOFTWARE/anaconda2/envs/vison/lib/python2.7/site-packages

dotask (*taskname, inputs*)

Generic test master function.

launchtask (*taskname*)

run (*explogf=None, elvis=None*)

```
wait_and_run (explogf, elvis='6.1.0')
```

3.1.2 FlatFielding.py

Flat-fielding Utilities.

Created on Fri Apr 22 16:13:22 2016

@author: raf

```
class vison.pipe.FlatFielding.FlatField (fitsfile='', data={}, meta={})
```

```
    parse_fits ()
```

```
vison.pipe.FlatFielding.fit2D (xx, yy, zz, degree=1)
```

```
vison.pipe.FlatFielding.get_ilum (img,    pdegree=5,    filtsize=15,    filtertype='median',  
                                Tests=False)
```

```
vison.pipe.FlatFielding.get_ilum_splines (img,    filtsize=25,    filtertype='median',  
                                         Tests=False)
```

```
vison.pipe.FlatFielding.produce_IndivFlats (infits, outfits, settings, runonTests, pro-  
                                         cesses=6)
```

```
vison.pipe.FlatFielding.produce_MasterFlat (infits, outfits, mask=None, settings={})
```

Produces a Master Flat out of a number of flat-illumination exposures. Takes the outputs from produce_IndivFlats.

```
vison.pipe.FlatFielding.produce_SingleFlatfield (infits, outfits, settings={}, runon-  
                                              Tests=False)
```

DATA MODEL

4.1 Data Model

4.1.1 ccd.py

Data model for Euclid-VIS CCDs (ground testing at MSSL)

Created on Fri Nov 13 17:42:36 2015

Author Ruyman Azzollini

class `vison.datamodel.ccd.CCD` (*infits=None, extensions=[-1], getallexensions=False*)

Class of CCD objects. Euclid Images as acquired by ELVIS software (Euclid LabView Imaging Software).

The class has been extended to handle multi-extension images. This is useful to also “host” calibration data-products, such as Flat-Fields.

add_extension (*data, header=None, label=None, headerdict=None*)

divide_by_flatfield (*FF, extension=-1*)

Divides by a Flat-field

get_cutout (*corners, Quadrant, canonical=False, extension=-1*)

Returns a cutout from the CCD image, either in canonical or non-canonical orientation.

Parameters

- **corners** (*list (of int)*) – [x0,x1,y0,y1]
- **Quadrant** (*char*) – Quadrant, one of ‘E’, ‘F’, ‘G’, ‘H’
- **canonical** (*bool*) – Canonical [True] = with readout-node at pixel index (0,0) regardless of quadrant. This is the orientation which corresponds to the data-readin order (useful for cross-talk measurements, for example). Non-Canonical [False] = with readout-node at corner matching placement of quadrant on the CCD. This is the orientation that would match the representation of the image on DS9.
- **extension** (*int*) – extension number. Default = -1 (last)

get_mask (*mask*)

get_quad (*Quadrant, canonical=False, extension=-1*)

Returns a quadrant in canonical or non-canonical orientation.

Parameters

- **Quadrant** (*char*) – Quadrant, one of ‘E’, ‘F’, ‘G’, ‘H’
- **canonical** –

Canonical [True] = with readout-node at pixel index (0,0) regardless of quadrant. This is the orientation which corresponds to the data-reading order (useful for cross-talk measurements, for example). Non-Canonical [False] = with readout-node at corner matching placement of quadrant on the CCD. This is the orientation that would match the representation of the image on DS9.

Parameters **extension** (*int*) – extension number. Default = -1 (last)

get_stats (*Quadrant*, *sector*='img', *statkeys*=['mean'], *trimscan*=[0, 0], *extension*=-1)

getsectioncollims (*QUAD*)

Returns limits of sections: prescan, image and overscan

set_quad (*inQdata*, *Quadrant*, *canonical*=False, *extension*=-1)

sim_window (*vstart*, *vend*, *extension*=-1)

simadd_flatilum (*levels*={'H': 0.0, 'E': 0.0, 'G': 0.0, 'F': 0.0}, *extension*=-1)

simadd_points (*flux*, *fwhm*, *CCDID*='CCD1', *dx*=0, *dy*=0, *extension*=-1)

simadd_poisson (*extension*=-1)

simadd_ron (*extension*=-1)

sub_bias (*superbias*, *extension*=-1)

Subtracts a superbias

sub_offset (*Quad*, *method*='row', *scan*='pre', *trimscan*=[3, 2], *extension*=-1)

writeto (*fitsf*, *clobber*=False, *unsigned16bit*=False)

class vison.datamodel.ccd.**Extension** (*data*, *header*=None, *label*=None, *headerdict*=None)

Extension Class

vison.datamodel.ccd.**test_create_from_scratch**()

vison.datamodel.ccd.**test_load_ELVIS_fits**()

4.1.2 EXPLOGtools.py

class vison.datamodel.EXPLOGtools.**ExpLogClass** (*elvis*='5.7.04')

addRow (*row*)

iniExplog ()

writeto (*outfile*)

vison.datamodel.EXPLOGtools.**iniExplog** (*elvis*)

vison.datamodel.EXPLOGtools.**loadExpLog** (*expfile*, *elvis*='5.7.04')

Loads an Exposure Log from file.

vison.datamodel.EXPLOGtools.**mergeExpLogs** (*explogList*, *addpedigree*=False)

Merges explog objects in a list.

vison.datamodel.EXPLOGtools.**test** ()

This Tests needs UPDATE (for data access and probably data format)

4.1.3 HKtools.py

House-Keeping inspection and handling tools.

History

Created on Thu Mar 10 12:11:58 2016

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

`vison.datamodel.HKtools.HKplot` (*allHKdata*, *keylist*, *key*, *dtobjs*, *filename*='', *stat*='mean')

Plots the values of a HK parameter as a function of time.

Parameters

- **allHKdata** – HKdata = [(nfiles,nstats,nHKparams)]
- **keylist** – list with all HK keys.
- **key** – selected key.
- **tdeltahour** – time axis.

Returns None!!

`vison.datamodel.HKtools.filtervalues` (*values*, *key*)

`vison.datamodel.HKtools.iniHK_QFM` (*elvis*='6.0.0', *length*=0)

`vison.datamodel.HKtools.loadHK_QFM` (*filename*, *elvis*='5.8.X')

Loads a HK file

Structure: tab separated columns, one per Keyword. First column is a timestamp, and there may be a variable number of rows (readings).

Parameters

- **filename** – path to the file to be loaded, including the file itself
- **form** – format of HK file, given by version of “ELVIS”

Returns dictionary with pairs parameter:[values]

`vison.datamodel.HKtools.loadHK_preQM` (*filename*, *elvis*='5.7.07')

Loads a HK file

It only assumes a structure given by a HK keyword followed by a number of of tab-separated values (number not specified). Note that the length of the values arrays is variable (depends on length of exposure and HK sampling rate).

Parameters

- **filename** – path to the file to be loaded, including the file itself
- **form** – format of HK file, given by version of “ELVIS”

Returns dictionary with pairs parameter:[values]

`vison.datamodel.HKtools.parseDTstr` (*DTstr*)

`vison.datamodel.HKtools.parseHKfiles` (*HKlist*, *elvis*='5.7.07')

Parameters **HKlist** – list of HK files (path+name).

Returns [obsids],[dtobjs],[tdeltasec],[HK_keys], [data(nfiles,nstats,nHKparams)]

`vison.datamodel.HKtools.parseHKfname (HKfname)`

Parses name of a HK file to retrieve OBSID, date and time, and ROE number.

Parameters `HKfname` – name of HK file.

Returns `obsid,dtobj=datetime.datetime(yy,MM,dd,hh,mm,ss),ROE`

`vison.datamodel.HKtools.reportHK (HKs, key, reqstat='all')`

Returns (mean, std, min, max) for each keyword in a list of HK dictionaries (output from loadHK).

Parameters

- **HK** – dictionary with HK data.
- **key** – HK key.

Reqstat what statistic to retrieve.

`vison.datamodel.HKtools.synthHK (HK)`

Synthesizes the values for each parameter in a HK dictionary into [mean,std,min,max].

Parameters `HK` – a dictionary as those output by loadHK.

Returns dictionary with pairs parameter:[mean,std,min,max]

4.1.4 QLAtools.py

Quick-Look-Analysis Tools.

History

Created on Wed Mar 16 11:31:58 2016

@author: Ruyman Azzollini

`vison.datamodel.QLAtools.dissectFITS (FITSfile, path='')`

`vison.datamodel.QLAtools.getacrosscolscut (CCDobj)`

`vison.datamodel.QLAtools.getacrossrowscut (CCDobj)`

`vison.datamodel.QLAtools.getsectionstats (CCDobj, QUAD, section, xbuffer=(0, 0),
ybuffer=(0, 0))`

`vison.datamodel.QLAtools.plotAcCOLcuts (dissection, filename=None, suptitle='')`

`vison.datamodel.QLAtools.plotAcROWcuts (dissection, filename=None, suptitle='')`

`vison.datamodel.QLAtools.plotQuads (CCDobj, filename=None, suptitle='')`

`vison.datamodel.QLAtools.reportFITS (FITSfile, outpath='')`

ANALYSIS (SHARED)

5.1 Analysis (Shared)

5.1.1 ellipse.py

Auxiliary module with functions to generate generalized ellipse masks.

author Ruyman Azzollini

contact r.azzollini@ucl.ac.uk

`vison.analysis.ellipse.area_superellip` (*r, q, c=0*)

Returns area of superellipse, given the semi-major axis length

`vison.analysis.ellipse.dist_superellipse` (*n, center, q=1, pos_ang=0.0, c=0.0*)

Form an array in which the value of each element is equal to the semi-major axis of the superellipse of specified center, axial ratio, position angle, and c parameter which passes through that element. Useful for super-elliptical aperture photometry.

Inspired on `dist_ellipse.pro` from AstroLib (IDL).

Note: this program doesn't take into account the change in the order of axes from IDL to Python. That means, that in 'n' and in 'center', the order of the coordinates must be reversed with respect to the case for `dist_ellipse.pro`, in order to get expected results. Nonetheless, the polar angle means the counter-clock wise angle with respect to the 'y' axis.

Parameters

- **n** – shape of array (N1,N2)
- **center** – center of superellipse radii: (c1,c2)
- **q** – axis ratio r2/r1
- **pos_ang** – position angle of isophotes, in degrees, CCW from axis 1
- **c** – boxyness (*c*>0) /diskyness (*c*<0)

`vison.analysis.ellipse.effective_radius` (*area, q=0, c=0*)

Returns semi-major axis length of superellipse, given the area

5.1.2 Guyonnet15.py

Library with functions that implement the algorithms described in Guyonnet+15. "Evidence for self-interaction of charge distribution in CCDs" Guyonnet, Astier, Antilogus, Regnault and Doherty 2015

Notes:

- I renamed “x” (pixel boundary index) to “b”, to avoid confusion with cartesian “x”.
- In paper, X belongsto [(0,1),(1,0),(0,-1),(-1,0)]. Here b is referred to as cardinal points “N”, “E”, “S”, “W”. It is linked to matrix index ib, running between 0 and 3.

Created on Thu Sep 22 11:38:24 2016

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

`vison.analysis.Guyonnet15.correct_estatic` (*img*, *aijb*)

Corrects an image from pixel-boundaries deformation due to electrostatic forces. Subtracts delta-Q.

Parameters

- **img** – image, 2D array
- **aijb** – Aijb matrix, 3D array

Returns array, img - delta-Q

`vison.analysis.Guyonnet15.degrade_estatic` (*img*, *aijb*)

Degrades an image according to matrix of pixel-boundaries deformations. Follows on Eq. 11 of G15. Adds delta-Q.

Parameters

- **img** – image, 2D array
- **aijb** – Aijb matrix, 3D array

Returns array, img + delta-Q

`vison.analysis.Guyonnet15.fpred_aijb` (*p*, *i*, *j*, *ib*)

‘The smoothing model assumes that a_{ij}^x coefficients are the product of a function of distance from the source charge to the considered boundary (r_{ij}) and that it also trivially depends on the angle between the source-boundary vector and the normal to the boundary (θ_{ij}^x)’

Eq. 18

Parameters

- **p** – parameters of the radial function (list of 2)
- **i** – pixel coordinate i
- **j** – pixel coordinate j
- **ib** – boundary index [0, 1, 2, 3]

Returns $f(r_{ij})\cos(\theta_{ij}^x)$

`vison.analysis.Guyonnet15.frdist` (*i*, *j*, *ib*)

Distance from the source charge to considered boundary “b”

Parameters

- **i** – pixel coordinate i
- **j** – pixel coordinate j
- **ib** – boundary index [0, 1, 2, 3]

Returns distance $r(ijb)$

`vison.analysis.Guyonnet15.ftheta_bound` (*i*, *j*, *ib*)

“ θ_{ij}^x is] the angle between the source-boundary vector and the normal to the boundary”.

Parameters

- **i** – pixel coordinate i
- **j** – pixel coordinate j
- **ib** – boundary index [0, 1, 2, 3]

Returns $\theta_{i,j}^x$

`vison.analysis.Guyonnet15.fun_p(x, *p)`
 auxiliary function to 'solve_for_psmooth'

`vison.analysis.Guyonnet15.generate_GaussPSF(N, sigma)`
 Create a circular symmetric Gaussian centered on the centre of a NxN matrix/image.

`vison.analysis.Guyonnet15.get_Rdisp(img, aijb)`
 Retrieves map of relative displacements of pixel boundaries, for input img and Aijb matrix.
 See G15 - Eq. 6

Parameters

- **img** – image, 2D array
- **aijb** – aijb matrix, 3D array NxNx4

Returns array, relative displacements all boundaries of pixels in img

`vison.analysis.Guyonnet15.get_deltaQ(img, aijb, writeFits=False)`
 Retrieves deltaQ map for input image and aijb matrix.

See G15 - Eq. 11

Parameters

- **img** – image, 2D array
- **aijb** – Aijb matrix, 3D array
- **writeFits** – save FITS file with resulting dQ map (optional)

Returns array, matrix with delta-Q for each pixel in img, given aijb

`vison.analysis.Guyonnet15.get_kernel(aijb, writeFits=False)`
 'kernel' is an array (2N-1)x(2N-1)x4. Each plane kernel[:,b] is a 2D array with the displacement coefficients aijb, in all directions around a pixel at (0,0).

Parameters

- **aijb** – array, matrix with displacements in 1st quadrant
- **writeFits** – save kernel to 4 FITS files

Returns kernel matrix, (2N-1)x(2N-1)x4

`vison.analysis.Guyonnet15.plot_map(z, ii, jj, title='')`

`vison.analysis.Guyonnet15.plot_maps_ftheta(f, ii, jj, suptitle='')`

`vison.analysis.Guyonnet15.show_disps_CCD273(aijb, stretch=5.0, peak=28571.428571428572, N=25, sigma=1.6, title='', filename='')`

```
vison.analysis.Guyonnet15.solve_for_A_linalg(covij, var=1.0, mu=1.0, doplot=False,
                                             psmooth=None, returnAll=False)
```

Function to retrieve the A matrix of pixel boundaries displacements, given a matrix of pixel covariances, variance, and mu.

if var==1 and mu==1, it is understood that covij is the correlation matrix.

See section 6.1 of G15.

Parameters

- **covij** – array, squared matrix with pixel covariances.
- **var** – float, variance of the flat-field.
- **mu** – float, mean value of the flat-field.
- **doplot** – if True, plot the fit of the fpred(ijb) function
- **psmooth** – coefficients of the fpred(aijb) function (Eq. 18)
- **returnAll** – bool, controls return values

Returns if returnAll == True, return (aijb, psmooth), otherwise return aijb only

```
vison.analysis.Guyonnet15.solve_for_psmooth(covij, var, mu, doplot=False)
```

Solving (p0,p1) parameters in Eq. 18 using covariance matrix and measured covariance matrix.

Parameters

- **covij** – array, covariance matrix
- **var** – float, variance
- **mu** – float, expected value of pixel values (“mean” of flat-field)
- **doplot** – bool, if True, plot data and best fit model

Returns best-fit parameters, and errors: 2 tuples of 2 elements each

```
vison.analysis.Guyonnet15.test0()
```

```
vison.analysis.Guyonnet15.test_getkernel()
```

```
vison.analysis.Guyonnet15.test_selfconsist()
```

```
vison.analysis.Guyonnet15.test_solve()
```

CHARGE INJECTION TOOLS

6.1 Charge Injection Tools

6.1.1 lib.py

NEEDSREVISION

Module to provide common tools for analysis of Charge Injection acquisitions.

Created on Thu Sep 14 15:32:10 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

6.1.2 plot.py

NEEDSREVISION

Charge Injection Plotting Tools.

Created on Thu Sep 14 15:39:34 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

“FLAT” ACQ. ANALYSIS TOOLS

7.1 “Flat” Acq. Analysis Tools

7.1.1 `ptc.py`

NEEDSREVISION

Module with tools used in PTC analysis.

Created on Thu Sep 14 16:29:36 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

`vison.flat.ptc.fitPTC` (*means, var*)

8.1 Image Analysis

8.1.1 bits.py

NEEDSREVISION

Image bits analysis tools.

Created on Thu Sep 14 15:54:14 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

MONITORING (“EYEGORE”)

9.1 Monitoring (“Eyegore”)



Fig. 9.1: You must be Igor...

9.1.1 eyegore.py

eyegore

data acquisition monitoring script for vison package.

‘- You must be Igor... - No, it’s pronounced “Eye-gore”.’

Created on Thu Feb 2 15:27:39 2017

Author Ruyman Azzollini

```
class vison.eyegore.eyegore.ExpLogDisplay (parent, path, elvis='6.0.0')
```

```
changeNumeric (data)
    if the data to be sorted is numeric change to float

get_data ()

isNumeric (s)
    test if a string s is numeric

search_EXPLOG ()

sortBy (tree, col, descending)
    sort tree contents when a column header is clicked

class vison.eyegore.eyegore.Eyegore (path, interval=3000)

class vison.eyegore.eyegore.HKDisplay (parent, path, elvis= '6.0.0')

    get_data ()

    search_HKfile ()

    select_HKkeys ()

class vison.eyegore.eyegore.ImageDisplay (parent, path)

    gen_render ()

    get_data ()
        PENDING
```

10.1 OGSE Tools

10.1.1 ogse.py

EUCLID-VIS Ground Calibration Campaign

Model of the calibration OGSE

Created on Fri Sep 8 12:11:55 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

`vison.ogse.ogse.get_FW_ID(wavelength)`

returns FW key corresponding to input wavelength. :param wavelength: integer, wavelength.

POINT-SOURCE ANALYSIS

11.1 Point-Source Analysis

11.1.1 basis.py

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

Created on Thu Apr 20 18:56:40 2017

```
class vison.point.basis.SpotBase (data, log=None)
```

11.1.2 display.py

Display Library for Point-Source Analysis

requires matplotlib

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

Created on Fri Apr 21 14:02:57 2017

```
vison.point.display.show_spots_allCCDs (spots_bag, title='', filename='', dobar=True)
```

11.1.3 gauss.py

Gaussian Model of Point-like Sources

Simple class to do Gaussian Fitting to a spot.

requires NumPy, astropy

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

Created on Thu Apr 20 16:42:47 2017

```
class vison.point.gauss.Gaussmeter (data, log=None, **kwargs)
```

Provides methods to measure the shape of an object using a 2D Gaussian Model.

Parameters

- **data** (*np.ndarray*) – stamp to be analysed.
- **log** (*instance*) – logger
- **kwargs** (*dict*) – additional keyword arguments

Settings dictionary contains all parameter values needed.

fit_Gauss ()

11.1.4 models.py

Models (Point-Like Sources)

Library module with models for processing of point-source imaging data.

requires NumPy

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

Created on Wed Apr 19 11:47:00 2017

`vison.point.models.fgauss2D(x, y, p)`

A gaussian fitting function where $p[0]$ = amplitude $p[1]$ = x_0 $p[2]$ = y_0 $p[3]$ = σ_x $p[4]$ = σ_y $p[5]$ = floor

11.1.5 photom.py

Aperture Photometry of point-like objects

Simple class to do aperture photometry on a stamp of a point-source.

requires NumPy

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

Created on Thu Apr 20 14:37:46 2017

class `vison.point.photom.Photometer` (*data*, *log=None*, ***kwargs*)

Provides methods to measure the shape of an object.

Parameters

- **data** (*np.ndarray*) – stamp to be analysed.
- **log** (*instance*) – logger
- **kwargs** (*dict*) – additional keyword arguments

Settings dictionary contains all parameter values needed.

doap_photom (*centre*, *rap*, *rin=-1.0*, *rout=-1.0*, *gain=3.5*, *doErrors=True*, *subbgd=False*)

get_centroid (*rap=None*, *full=False*)

TODO: add aperture masking

measure_bgd (*rin*, *rout*)

sub_bgd (*rin*, *rout*)

11.1.6 shape.py

Quadrupole Moments Shape Measurement

Simple class to measure quadrupole moments and ellipticity of an object.

requires NumPy, PyFITS

author Sami-Matias Niemi, Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

class `vison.point.shape.Shapemeter` (*data*, *log=None*, ***kwargs*)
Provides methods to measure the shape of an object.

Parameters

- **data** (*np.ndarray*) – stamp to be analysed.
- **log** (*instance*) – logger
- **kwargs** (*dict*) – additional keyword arguments

Settings dictionary contains all parameter values needed.

circular2DGaussian (*x*, *y*, *sigma*)
Create a circular symmetric Gaussian centered on x, y.

Parameters

- **x** (*float*) – x coordinate of the centre
- **y** (*float*) – y coordinate of the centre
- **sigma** (*float*) – standard deviation of the Gaussian, note that `sigma_x = sigma_y = sigma`

Returns circular Gaussian 2D profile and x and y mesh grid

Return type dict

ellip2DGaussian (*x*, *y*, *sigmax*, *sigmay*)
Create a two-dimensional Gaussian centered on x, y.

Parameters

- **x** (*float*) – x coordinate of the centre
- **y** (*float*) – y coordinate of the centre
- **sigmax** (*float*) – standard deviation of the Gaussian in x-direction
- **sigmay** (*float*) – standard deviation of the Gaussian in y-direction

Returns circular Gaussian 2D profile and x and y mesh grid

Return type dict

measureRefinedEllipticity ()
Derive a refined iterated polarisability/ellipticity measurement for a given object.

By default polarisability/ellipticity is defined in terms of the Gaussian weighted quadrupole moments. If `self.shsettings['weighted']` is False then no weighting scheme is used.

The number of iterations is defined in `self.shsettings['iterations']`.

Returns centroids [indexing stars from 1], ellipticity (including projected e1 and e2), and R2

Return type dict

quadrupoles (*image*)

Derive quadrupole moments and ellipticity from the input image.

Parameters **img** (*ndarray*) – input image data

Returns quadrupoles, centroid, and ellipticity (also the projected components e1, e2)

Return type dict

writeFITS (*data*, *output*)

Write out a FITS file using PyFITS.

Parameters

- **data** (*ndarray*) – data to write to a FITS file
- **output** (*string*) – name of the output file

Returns None

11.1.7 spot.py

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

Created on Thu Apr 20 15:35:08 2017

class `vison.point.spot.Spot` (*data*, *log=None*, ***kwargs*)

Provides methods to do point-source analysis on a stamp. Aimed at basic analysis:

- Photometry
- Quadrupole Moments
- Gaussian Fit

Parameters

- **data** (*np.ndarray*) – stamp to be analysed.
- **log** (*instance*) – logger
- **kwargs** (*dict*) – additional keyword arguments

Settings dictionary contains all parameter values needed.

11.1.8 lib.py

FM-Calib. Campaign.

Library module with useful data and functions for processing of point-source imaging data.

Created on Wed Apr 5 10:21:05 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

```
vison.point.lib.fwcentroid(image, checkbox=1, maxiterations=30, threshold=1e-05,  
                           halfwidth=35, verbose=False, full=False)
```

Implement the Floating-window first moment centroid algorithm chosen for JWST target acquisition.

See JWST-STScI-001117 and JWST-STScI-001134 for details.

This code makes no attempt to vectorize or optimize for speed; it's pretty much just a straight verbatim implementation of the IDL-like pseudocode provided in JWST-STScI-001117

image [array_like] image to centroid

checkbox [int] size of moving checkbox for initial peak pixel guess. Default 1

maxiterations [int] Max number of loops. Default 30

halfwidth [int] Half width of the centroid box size (less 1). Specify as a scalar, or a tuple Xhalfwidth, Yhalfwidth. Empirical tests suggest this parameter should be at *least* the PSF FWHM for convergence, preferably some small factor larger

threshold [float] Position threshold for convergence

(ycen, xcen) [float tuple] Measured centroid position. Note that this is returned in Pythonic Y,X order for use as array indices, etc.

-Marshall Perrin 2011-02-11

```
vison.point.lib.get_Gauss_spot()
```

measurements: 'i0_G', 'ei0_G', 'flu_G', 'eflu_G', 'x0_G', 'ex0_G', 'y0_G', 'ey0_G',
'sx_G', 'esx_G', 'sy_G', 'esy_G', 'e1', 'e2', 'e', 'R2'

```
vison.point.lib.get_Quadrupoles_spot()
```

measurements: 'e1', 'e2', 'e', 'R2'

```
vison.point.lib.get_photom_spot(stamp, x0, y0)
```

measurements: 'apflu', 'eapflu', 'bgd', 'ebgd'

TEST SCRIPTS

12.1 Charge Injection Scripts

12.1.1 Charge Injection Scripts

CHINJ01

VIS Ground Calibration TEST: CHINJ01

Charge injection calibration (part 1) Injection vs. IG1-IG2

Created on Tue Aug 29 17:36:00 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

`vison.inject.CHINJ01.build_CHINJ01_scriptdict` (*IDL, IDH, IG1s, id_delays, toi_chinj, diff-values={}, elvis='6.0.0'*)

Builds CHINJ01 script structure dictionary.

Parameters

- **IDL** – int, [mV], value of IDL (Inject. Drain Low).
- **IDH** – int, [mV], Injection Drain High.
- **IG1s** – list of 2 ints, [mV], [min,max] values of IG1.
- **id_delays** – list of 2 ints, [mV], injection drain delays (2).
- **toi_chinj** – int, [us], TOI-charge injection.
- **diffvalues** – dict, opt, differential values.

CHINJ02

VIS Ground Calibration TEST: CHINJ02

Charge injection calibration (part 2) Injection vs. IDL (injection threshold)

Created on Tue Aug 29 17:36:00 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

`vison.inject.CHINJ02.build_CHINJ02_scriptdict` (*IDLs, IDH, id_delays, toi_chinj, diffvalues={}, elvis='6.0.0'*)

Builds CHINJ02 script structure dictionary.

Parameters

- **IDLs** – list of 2 ints, [mV], [min,max] values of IDL (Inject. Drain Low).
- **IDH** – int, [mV], Injection Drain High.
- **id_delays** – list of 2 ints, [mV], injection drain delays (2).
- **toi_chinj** – int, [us], TOI-charge injection.
- **diffvalues** – dict, opt, differential values.

12.2 Dark Scripts

12.2.1 “Dark Acquisitions” Scripts

BIAS01

VIS Ground Calibration TEST: BIAS01

Bias-structure/RON analysis script

Created on Tue Aug 29 16:53:40 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

`vison.dark.BIAS01.basic_analysis` (*DataDict, report, inputs, log=None*)
METACODE

Basic analysis of data.

`vison.dark.BIAS01.build_BIAS01_scriptdict` (*N, diffvalues={}, elvis='6.0.0'*)
Builds BIAS01 script structure dictionary.

Parameters

- **N** – integer, number of frames to acquire.
- **diffvalues** – dict, opt, differential values.
- **elvis** – char, ELVIS version.

`vison.dark.BIAS01.check_data` (*DataDict, report, inputs, log=None*)
METACODE

Checks quality of ingested data.

TODO: consider to raise an exception that would halt execution of task if processing data could be just a waste of time.

TODO: consider add a common binary “flags” variable as input/output. It could go in DataDict, and reported in log and report.

`vison.dark.BIAS01.feeder` (*inputs, elvis='6.1.0'*)

`vison.dark.BIAS01.filterexposures` (*structure, explogf, datapath, OBSID_lims, elvis*)

`vison.dark.BIAS01.meta_analysis` (*DataDict, report, inputs, log=None*)

`vison.dark.BIAS01.prep_data` (*DataDict*, *report*, *inputs*, *log=None*)
METACODE

Preparation of data for further analysis.

DARK01

VIS Ground Calibration TEST: DARK01

“Dark Current” analysis script

Created on Tue Aug 29 17:21:00 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

`vison.dark.DARK01.basic_analysis` (*DataDict*, *report*, *inputs*, *log=None*)
METACODE

Basic analysis of data.

`vison.dark.DARK01.build_DARK01_scriptdict` (*N*, *exptime*, *diffvalues={}*, *elvis='6.0.0'*)
Builds DARK01 script structure dictionary.

Parameters

- **N** – integer, number of frames to acquire.
- **exptime** – integer, ms, exposure time.
- **diffvalues** – dict, opt, differential values.

`vison.dark.DARK01.check_data` (*DataDict*, *report*, *inputs*, *log=None*)
METACODE

Checks quality of ingested data.

`vison.dark.DARK01.meta_analysis` (*DataDict*, *report*, *inputs*, *log=None*)

`vison.dark.DARK01.prep_data` (*DataDict*, *report*, *inputs*, *log=None*)
METACODE

Preparation of data for further analysis.

12.3 Flat-Illumination Scripts

12.3.1 Flat-Illumination Scripts

FLAT0X

VIS Ground Calibration TEST: FLAT0X

Flat-fields acquisition / analysis script

Created on Tue Aug 29 17:32:52 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

```
vison.flat.FLAT0X.build_FLAT0X_scriptdict (exptimes, wavelength=800, testkey='FLAT0X',  
                                             diffvalues={}, elvis='6.0.0')
```

Builds FLAT0X script structure dictionary.

Parameters

- **exptimes** – list of ints, exposure times.
- **wavelength** – int, wavelength.
- **testkey** – char, test identifier.
- **diffvalues** – dict, opt, differential values.

NL01

VIS Ground Calibration TEST: NL01

End-To-End Non-Linearity Curve

Tasks:

- Select exposures, get file names, get metadata (commandig, HK).
- Check exposure time pattern matches test design.
- Check quality of data (rough scaling of fluences with Exposure times).
- Subtract offset level.
- Divide by Flat-field.
- **Synoptic analysis:** fluence ratios vs. exptime ratios >> non-linearity curve
- extract: Non-Linearity curve for each CCD and quadrant
- produce synoptic figures
- Save results.

Created on Mon Apr 3 17:38:00 2017

author raf

contact r.azzollini_at_ucl.ac.uk

```
vison.flat.NL01.basic_analysis_NL01 (DataDict, RepDict, inputs, log=None)
```

Performs basic analysis: - builds Non-Lin curves - fits Non-Linearity curves

```
vison.flat.NL01.build_NL01_scriptdict (expts, exptinter, frames, wavelength=0, diffvalues={},  
                                         elvis='6.0.0')
```

Builds NL01 script structure dictionary.

Parameters

- **expts** – list of ints [ms], exposure times.
- **exptinter** – int, ms, exposure time of interleaved source-stability exposures.
- **frames** – list of ints, number of frames for each exposure time.
- **wavelength** – int, wavelength. Default: 0 (Neutral Density Filter)
- **diffvalues** – dict, opt, differential values.

```
vison.flat.NL01.filterexposures_NLC01(inwavelength, explogf, datapath, OBSID_lims, structure={ 'Ncols': 6, 'col6': { 'Exptime': 18.0, 'N': 2}, 'col4': { 'Exptime': 10.0, 'N': 2}, 'col5': { 'Exptime': 15.0, 'N': 2}, 'col2': { 'Exptime': 1.0, 'N': 2}, 'col3': { 'Exptime': 5.0, 'N': 2}, 'col1': { 'Exptime': 0, 'N': 5}}, elvis='5.7.04'))
```

Loads a list of Exposure Logs and selects exposures from test PSF0X.

The filtering takes into account an expected structure for the acquisition script.

The datapath becomes another column in DataDict. This helps dealing with tests that run overnight and for which the input data is in several date-folders.

```
vison.flat.NL01.prep_data_NL01(DataDict, RepDict, inputs, log=None)
```

Takes Raw Data and prepares it for further analysis. Also checks that data quality is enough.

```
vison.flat.NL01.run(inputs, log=None)
```

Test NL01 master function.

PTC0X

VIS Ground Calibration TEST: PTC_0X

Photon-Transfer-Curve Analysis PTC01 - nominal temperature PTC02 - alternative temperatures

Tasks:

- Select exposures, get file names, get metadata (commandig, HK).
- Check exposure time pattern matches test design.
- Check quality of data (rough scaling of fluences with Exposure times).
- Subtract pairs of exposures with equal fluence
- **Synoptic analysis:** variance vs. fluence variance(binned difference-frames) vs. fluence
- extract: RON, gain, gain(fluence)
- produce synoptic figures
- Save results.

Created on Mon Apr 3 17:00:24 2017

author raf

contact r.azzollini_at_ucl.ac.uk

```
vison.flat.PTC0X.basic_analysis_PTC0X(DataDict, RepDict, inputs, log=None)
```

Performs basic analysis of images: - builds PTC curves: both on non-binned and binned images

```
vison.flat.PTC0X.build_PTC0X_scriptdict(exptimes, frames, wavelength=800, diffvalues={}, elvis='6.0.0')
```

Builds PTC0X script structure dictionary.

Parameters

- **exptimes** – list of ints [ms], exposure times.
- **frames** – list of ints, number of frames for each exposure time.
- **wavelength** – int, wavelength. Default: 800 nm.
- **diffvalues** – dict, opt, differential values.

```
vison.flat.PTC0X.filterexposures_PTC0X(inwavelength, explogf, datapath, OBSID_lims, structure={ 'Ncols': 6, 'col6': { 'Exptime': 18.0, 'N': 2}, 'col4': { 'Exptime': 10.0, 'N': 2}, 'col5': { 'Exptime': 15.0, 'N': 2}, 'col2': { 'Exptime': 1.0, 'N': 2}, 'col3': { 'Exptime': 5.0, 'N': 2}, 'col1': { 'Exptime': 0, 'N': 5}}, elvis='5.7.04'))
```

Loads a list of Exposure Logs and selects exposures from test PSF0X.

The filtering takes into account an expected structure for the acquisition script.

The datapath becomes another column in DataDict. This helps dealing with tests that run overnight and for which the input data is in several date-folders.

```
vison.flat.PTC0X.meta_analysis_PSF0X(DataDict, RepDict, inputs, log=None)
```

Analyzes the variance and fluence: gain, and gain(fluence)

```
vison.flat.PTC0X.prep_data_PTC0X(DataDict, RepDict, inputs, log=None)
```

Takes Raw Data and prepares it for further analysis. Also checks that data quality is enough.

```
vison.flat.PTC0X.run(inputs, log=None)
```

Test PTC0X master function.

12.4 Other Test Scripts

12.4.1 Other Scripts

PERSIST01

VIS Ground Calibration TEST: PERSIST01

CCD Persistence test

Created on Tue Aug 29 17:39:00 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

```
vison.other.PERSIST01.build_PER01_scriptdict(exptSATUR, exptLATEN, diffvalues={}, elvis='6.0.0')
```

Builds PERSISTENCE01 script structure dictionary.

Parameters

- **exptSATUR** – int, saturation exposure time.
- **exptLATEN** – int, latency exposure time.
- **diffvalues** – dict, opt, differential values.

12.5 Point-Source Scripts

12.5.1 Point-Source Scripts

FOCUS00

VIS Ground Calibration TEST: FOCUS00

Focus analysis script

Tasks:

- Select exposures, get file names, get metadata (commandig, HK).
- Check quality of data (integrated fluxes are roughly constant, matching expected level).
- Subtract offset level.
- Divide by Flat-field.
- **Crop stamps of the sources on each CCD/Quadrant.**
 - save snapshot figures of sources.
- **for each source (5 x Nquadrants):**
 - measure shape using Gaussian Fit
- Find position of mirror that minimizes PSF sizes
- **Produce synoptic figures:** source size and ellipticity across combined FOV (of 3 CCDs)
- Save results.

Created on Mon Apr 03 16:21:00 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

```
vison.point.FOCUS00.basic_analysis_FOCUS00 (DataDict, Report, inputs, log=None, debug=False)
```

Performs basic analysis on spots: - 2D Gaussian Model shape measurements - Quadrupole Moments shape measurements

```
vison.point.FOCUS00.build_FOCUS00_scriptdict (wavelength, exptime, diffvalues={}, elvis='6.0.0')
```

Builds FOCUS00 script structure dictionary.

Parameters

- **wavelength** – int, [nm], wavelength.
- **exptime** – int, [ms], exposure time.
- **diffvalues** – dict, opt, differential values.

```
vison.point.FOCUS00.filterexposures_FOCUS00 (inwavelength, explogf, datapath, OB-SID_lims, structure={'Ncols': 6, 'col6': {'Mirr_pos': 70.2, 'Exptime': 10.0, 'N': 2}, 'col4': {'Mirr_pos': 70.0, 'Exptime': 10.0, 'N': 2}, 'col5': {'Mirr_pos': 70.1, 'Exptime': 10.0, 'N': 2}, 'col2': {'Mirr_pos': 69.8, 'Exptime': 10.0, 'N': 2}, 'col3': {'Mirr_pos': 69.9, 'Exptime': 10.0, 'N': 2}, 'col1': {'Mirr_pos': 69.8, 'Exptime': 0, 'N': 5}}, elvis='5.7.04')
```

Loads a list of Exposure Logs and selects exposures from test FOCUS00.

The filtering takes into account an expected structure for the acquisition script.

The datapath becomes another column in DataDict. This helps dealing with tests that run overnight and for which the input data is in several date-folders.

```
vison.point.FOCUS00.generate_Explog_FOCUS00 (wavelength, struct, elvis='6.0.0',
                                              date=datetime.datetime(1980, 2, 21, 7,
0))

vison.point.FOCUS00.generate_FITS_FOCUS00 (wavelength, explog, datapath, elvis='6.0.0')

vison.point.FOCUS00.generate_Fake_FOCUS00 (wavelength, date=datetime.datetime(1980, 2,
21, 7, 0), rootpath='')

    Generates Fake FOCUS00 data

vison.point.FOCUS00.generate_HK_FOCUS00 (explog, datapath, elvis='6.0.0')

vison.point.FOCUS00.get_FOCUS00_structure (wavelength)

vison.point.FOCUS00.get_basic_spot_FOCUS00 (stamp, x0, y0, log=None, debug=False)
    # TODO: # get basic statistics, measure and subtract background # update centroid # do aperture photometry #
    pack-up results and return

vison.point.FOCUS00.get_shape_spot_FOCUS00 (stamp, x0, y0, method='G', log=None, de-
bug=False)
    # TODO: # get basic statistics, measure and subtract background # update centroid # do aperture photometry #
    pack-up results and return

vison.point.FOCUS00.meta_analysis_FOCUS00 (DataDict, Report, inputs, log=None)
    Analyzes the relation between PSF shape and mirror position.
```

Parameters

- **DataDict** – Dictionary with input data
- **Report** – Report Objects
- **inputs** – Dictionary with inputs

```
vison.point.FOCUS00.prep_data_FOCUS00 (DataDict, Report, inputs, log=None, debug=False)
    Takes Raw Data and prepares it for further analysis. Also checks that data quality is enough.

vison.point.FOCUS00.run (inputs, log=None)
    Test FOCUS00 master function.
```

PSF0X

VIS Ground Calibration TEST: PSF0X

PSF vs. Fluence, and Wavelength PSF01 - nominal temperature PSF02 - alternative temperatures

Tasks:

- Select exposures, get file names, get metadata (commandig, HK).
- Check exposure time pattern matches test design.
- Check quality of data (rough scaling of fluences with Exposure times).
- Subtract offset level.
- Divide by Flat-field.
- **Crop stamps of the sources on each CCD/Quadrant.**
 - save snapshot figures of sources.
- **for each source:**
 - measure shape using weighted moments

- measure shape using Gaussian Fit
- Bayesian Forward Modelling the optomechanic+detector PSF
- Produce synoptic figures.
- Save results.

Created on Thu Dec 29 15:01:07 2016

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

`vison.point.PSF0X.basic_analysis_PSF0X` (*DataDict*, *RepDict*, *inputs*, *log=None*)

Performs basic analysis on spots: - Gaussian Fits: peak, position, width_x, width_y

`vison.point.PSF0X.bayes_analysis_PSF0X` (*DataDict*, *RepDict*, *inputs*, *log=None*)

Performs bayesian decomposition of the spot images:

- optomechanic PSF and detector PSF.

Also measures R2, fwhm and ellipticity of “extracted” detector PSF.

`vison.point.PSF0X.build_PSF0X_scriptdict` (*exptimes*, *frames*, *wavelength=800*, *diffvalues={}*, *elvis='6.0.0'*)

Builds PSF0X script structure dictionary.

Parameters

- **exptimes** – list of ints, [ms], exposure times.
- **frames** – list of frame numbers. Same length as exptimes.
- **wavelength** – int, [nm], wavelength.
- **diffvalues** – dict, opt, differential values.

`vison.point.PSF0X.filterexposures_PSF0X` (*inwavelength*, *explogf*, *datapath*, *OBSID_lims*, *structure={‘Ncols’: 6, ‘col6’: {‘Exptime’: 18.0, ‘N’: 3}, ‘col4’: {‘Exptime’: 10.0, ‘N’: 10}, ‘col5’: {‘Exptime’: 15.0, ‘N’: 4}, ‘col2’: {‘Exptime’: 1.0, ‘N’: 20}, ‘col3’: {‘Exptime’: 5.0, ‘N’: 18}, ‘col1’: {‘Exptime’: 0, ‘N’: 5}}, *elvis='5.7.04'*)*

Loads a list of Exposure Logs and selects exposures from test PSF0X.

The filtering takes into account an expected structure for the acquisition script.

The datapath becomes another column in DataDict. This helps dealing with tests that run overnight and for which the input data is in several date-folders.

`vison.point.PSF0X.meta_analysis_PSF0X` (*DataDict*, *RepDict*, *inputs*, *log=None*)

Analyzes the relation between detector PSF and fluence.

`vison.point.PSF0X.prep_data_PSF0X` (*DataDict*, *RepDict*, *inputs*, *log=None*)

Takes Raw Data and prepares it for further analysis. Also checks that data quality is enough.

`vison.point.PSF0X.run` (*inputs*, *log=None*)

Test PSF0X master function.

12.6 Trap-Pumping Scripts

12.6.1 Trap-Pumping Scripts

TP01

VIS Ground Calibration TEST: TP01

Trap-Pumping calibration (vertical)

Created on Tue Aug 29 17:37:00 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

```
vison.pump.TP01.build_TP01_scriptdict (Nshuffles_V, TOI_TPv, id_delays, diffvalues={},  
                                         elvis='6.0.0')
```

TP02

VIS Ground Calibration TEST: TP01

Trap-Pumping calibration (serial)

Created on Tue Aug 29 17:38:00 2017

author Ruyman Azzollini

contact r.azzollini__at__ucl.ac.uk

```
vison.pump.TP02.build_TP02_scriptdict (Nshuffles_H, dwell_hv, id_delays, diffvalues={},  
                                         elvis='6.0.0')
```

MISSING: different starting points (not implemented in ELVIS yet).

13.1 Scripts

13.1.1 HKmonitor.py

TODO: DEBUG, calls unexistent class LaTeX

Script to produce HK reports out of HK files in a folder. Aimed at quick inspection of data from Characterization and Calibration Campaigns of Euclid-VIS.

History

Created on Tue Mar 15 10:35:43 2016

@author: Ruyman Azzollini (MSSL)

13.1.2 quickds9.py

Wrap-up of ds9 to quickly load a number of images, for inspection.

History

Created on Thu Mar 17 13:18:10 2016

@author: Ruyman Azzollini

13.1.3 vis_genDataSet.py

EUCLID-VIS Ground Calibration Campaign

Development: Creating Calibration Campaign Fake Data-set

Created on Tue Sep 05 16:07:00 2017

autor Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

```
vison.scripts.vis_genDataSet.datasetGenerator (TestsSelector, doGenExplog, doGenHK,  
                                                doGenFITS, outpath, elvis, Nrows=0)
```

```
vison.scripts.vis_genDataSet.genExpLog (toGen, explogf, equipment, elvis='6.1.0')
```

13.1.4 vis_mkscripts.py

VIS Ground Calibration Campaign

Automatically Generating Calibration Campaign Scripts.

Created on Fri Sep 08 12:03:00 2017

autor Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

SUPPORT CODE

14.1 Support Code

NEEDS REVISION TODO: UPDATE

Module with functions related to the handling of ds9 from python through XPA.

History

Created on Thu Mar 17 13:18:10 2016

@author: Ruyman Azzollini

class `vison.support.ds9.ds9class`

A very simple class to handle ds9 through xpa.

isOpen ()

Return True if this ds9 window is open and available for communication, False otherwise.

launch ()

Launches ds9

xpaget (*cmd*)

Executes xpaget and retrieves the stdout. If an error happens, an exception is raised.

xpaset (*cmd*)

Executes xpaset.

zoomhere (*x, y, zoom*)

Zooms in on given coordinates of display (ds9).

Module to issue WARNING / ALERT phone calls to designated phone numbers.

‘... E.T. phone home...’

Created on Thu Sep 14 10:13:12 2017

@author: raf

class `vison.support.ET.ET`

Class to do phone calls.

dial_numbers (*url*)

Dials one or more phone numbers from a Twilio phone number.

Parameters *url* – char, URL with the TwiML code that Twilio uses as instructions on call.

Basically, it provides a message to be voiced, as intended.

`vison.support.ET.grab_numbers_and_codes` ()

Retrieves phone numbers and access codes necessary to make the phone calls.

IO related functions.

requires PyFITS

requires NumPy

author Sami-Matias Niemi

contact r.azzollini_at_ucl.ac.uk

`vison.support.files.cPickleDump` (*data*, *output*)
 Dumps data to a cPickled file.

Parameters

- **data** – a Python data container
- **output** – name of the output file

Returns None

`vison.support.files.cPickleDumpDictionary` (*dictionary*, *output*)
 Dumps a dictionary of data to a cPickled file.

Parameters

- **dictionary** – a Python data container does not have to be a dictionary
- **output** – name of the output file

Returns None

`vison.support.files.cPickleRead` (*file*)
 Loads data from a pickled file.

Euclid-VIS Calibration Programme Pipeline: vison

Reporting Utilities.

History

Created on Wed Jan 25 16:58:33 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

class `vison.support.report.Content` (*contenttype=''*)

class `vison.support.report.Figure` (*figpath*, *textfraction=0.7*, *caption=None*, *label=None*)

generate_Latex ()

Generates LaTeX as list of strings.

class `vison.support.report.Section` (*Title=''*, *level=0*)

generate_Latex ()

class `vison.support.report.Table` (*tableDict*, *formats={}*, *names=[]*, *caption=None*)

PENDING:

- adjust width of table to texwidth:

```
esizebox{ extwidth}{!}{
... end{tabular}}
```

- include option to rotate table to show in landscape

generate_Latex()

Generates LaTeX as list of strings.

class vison.support.report.**Text**(*text*)

generate_Latex()

Just a collection of LaTeX templates for use in report.py

History

Created on Mon Jan 30 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

vison.support.latex.**generate_header**(*test, model, author*)

These functions can be used for logging information.

Warning: logger is not multiprocessing safe.

author Sami-Matias Niemi

contact r.azzollini_at_ucl.ac.uk

version 0.3

class vison.support.logger.**SimpleLogger**(*filename, verbose=False*)

A simple class to create a log file or print the information on screen.

write(*text*)

Writes text either to file or screen.

vison.support.logger.**setUpLogger**(*log_filename, loggername='logger'*)

Sets up a logger.

Param log_filename: name of the file to save the log.

Param loggername: name of the logger

Returns logger instance

Euclid VIS Ground Calibration Campaign Classes and functions to do Real-Time Monitoring of the Data Acquisition

Created on Wed Feb 1 17:37:32 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

vison.support.monitor.**test**()

Tests Module Basic Functionality

Euclid-VIS Calibration Programme Pipeline: vison

Reporting Utilities.

History

Created on Wed Jan 25 16:58:33 2017

author Ruyman Azzollini

contact r.azzollini_at_ucl.ac.uk

class vison.support.report.**Content** (*contenttype=''*)

class vison.support.report.**Figure** (*figpath, textfraction=0.7, caption=None, label=None*)

generate_Latex ()

Generates LaTeX as list of strings.

class vison.support.report.**Section** (*Title='', level=0*)

generate_Latex ()

class vison.support.report.**Table** (*tableDict, formats={}, names=[], caption=None*)

PENDING:

- adjust width of table to texwidth:

esizebox{ *extwidth* }{!}{

... end{tabular}}

- include option to rotate table to show in landscape

generate_Latex ()

Generates LaTeX as list of strings.

class vison.support.report.**Text** (*text*)

generate_Latex ()

INDICES AND TABLES

- genindex
- modindex
- search

PYTHON MODULE INDEX

V

vison.analysis.ellipse, 15
vison.analysis.Guyonnet15, 15
vison.dark.BIAS01, 36
vison.dark.DARK01, 37
vison.datamodel.ccd, 11
vison.datamodel.EXPLOGtools, 12
vison.datamodel.HKtools, 13
vison.datamodel.QLAtools, 14
vison.eyegore.eyegore, 25
vison.flat.FLAT0X, 37
vison.flat.NL01, 38
vison.flat.ptc, 21
vison.flat.PTC0X, 39
vison.image.bits, 23
vison.inject.CHINJ01, 35
vison.inject.CHINJ02, 35
vison.inject.lib, 19
vison.inject.plot, 19
vison.ogse.ogse, 27
vison.other.PERSIST01, 40
vison.pipe.FlatFielding, 10
vison.pipe.master, 9
vison.point.basis, 29
vison.point.display, 29
vison.point.FOCUS00, 40
vison.point.gauss, 29
vison.point.lib, 32
vison.point.models, 30
vison.point.photom, 30
vison.point.PSF0X, 42
vison.point.shape, 31
vison.point.spot, 32
vison.pump.TP01, 44
vison.pump.TP02, 44
vison.scripts.HKmonitor, 45
vison.scripts.quickds9, 45
vison.scripts.vis_genDataSet, 45
vison.scripts.vis_mkscripts, 46
vison.support.ds9, 47
vison.support.ET, 47
vison.support.files, 48
vison.support.latex, 49
vison.support.logger, 49
vison.support.monitor, 49
vison.support.report, 49

A

add_extension() (vison.datamodel.ccd.CCD method), 11
 addRow() (vison.datamodel.EXPLOGtools.ExpLogClass method), 12
 area_superellip() (in module vison.analysis.ellipse), 15

B

basic_analysis() (in module vison.dark.BIAS01), 36
 basic_analysis() (in module vison.dark.DARK01), 37
 basic_analysis_FOCUS00() (in module vison.point.FOCUS00), 41
 basic_analysis_NL01() (in module vison.flat.NL01), 38
 basic_analysis_PSF0X() (in module vison.point.PSF0X), 43
 basic_analysis_PTC0X() (in module vison.flat.PTC0X), 39
 bayes_analysis_PSF0X() (in module vison.point.PSF0X), 43
 BIAS01 (vison.pipe.master.Pipe attribute), 9
 build_BIAS01_scriptdict() (in module vison.dark.BIAS01), 36
 build_CHINJ01_scriptdict() (in module vison.inject.CHINJ01), 35
 build_CHINJ02_scriptdict() (in module vison.inject.CHINJ02), 35
 build_DARK01_scriptdict() (in module vison.dark.DARK01), 37
 build_FLAT0X_scriptdict() (in module vison.flat.FLAT0X), 37
 build_FOCUS00_scriptdict() (in module vison.point.FOCUS00), 41
 build_NL01_scriptdict() (in module vison.flat.NL01), 38
 build_PER01_scriptdict() (in module vison.other.PERSIST01), 40
 build_PSF0X_scriptdict() (in module vison.point.PSF0X), 43
 build_PTC0X_scriptdict() (in module vison.flat.PTC0X), 39
 build_TP01_scriptdict() (in module vison.pump.TP01), 44
 build_TP02_scriptdict() (in module vison.pump.TP02), 44

C

CCD (class in vison.datamodel.ccd), 11
 changeNumeric() (vison.eyegore.eyegore.ExpLogDisplay method), 25
 check_data() (in module vison.dark.BIAS01), 36
 check_data() (in module vison.dark.DARK01), 37
 circular2DGaussian() (vison.point.shape.Shapemeter method), 31
 Content (class in vison.support.report), 48, 50
 correct_estatic() (in module vison.analysis.Guyonnet15), 16
 cPickleDump() (in module vison.support.files), 48
 cPickleDumpDictionary() (in module vison.support.files), 48
 cPickleRead() (in module vison.support.files), 48

D

datasetGenerator() (in module vison.scripts.vis_genDataSet), 45
 degrade_estatic() (in module vison.analysis.Guyonnet15), 16
 dial_numbers() (vison.support.ET.ET method), 47
 dissectFITS() (in module vison.datamodel.QLAtools), 14
 dist_superellipse() (in module vison.analysis.ellipse), 15
 divide_by_flatfield() (vison.datamodel.ccd.CCD method), 11
 doap_photom() (vison.point.photom.Photometer method), 30
 dotask() (vison.pipe.master.Pipe method), 9
 ds9class (class in vison.support.ds9), 47

E

effective_radius() (in module vison.analysis.ellipse), 15
 ellip2DGaussian() (vison.point.shape.Shapemeter method), 31
 ET (class in vison.support.ET), 47
 ExpLogClass (class in vison.datamodel.EXPLOGtools), 12
 ExpLogDisplay (class in vison.eyegore.eyegore), 25
 Extension (class in vison.datamodel.ccd), 12
 Eyegore (class in vison.eyegore.eyegore), 26

F

feeder() (in module vison.dark.BIAS01), 36
fgauss2D() (in module vison.point.models), 30
Figure (class in vison.support.report), 48, 50
filterexposures() (in module vison.dark.BIAS01), 36
filterexposures_FOCUS00() (in module vison.point.FOCUS00), 41
filterexposures_NLC01() (in module vison.flat.NL01), 38
filterexposures_PSF0X() (in module vison.point.PSF0X), 43
filterexposures_PTC0X() (in module vison.flat.PTC0X), 39
filtervalues() (in module vison.datamodel.HKtools), 13
fit2D() (in module vison.pipe.FlatFielding), 10
fit_Gauss() (vison.point.gauss.Gaussmeter method), 30
fitPTC() (in module vison.flat.ptc), 21
FlatField (class in vison.pipe.FlatFielding), 10
fpred_aijb() (in module vison.analysis.Guyonnet15), 16
frdlist() (in module vison.analysis.Guyonnet15), 16
ftheta_bound() (in module vison.analysis.Guyonnet15), 16
fun_p() (in module vison.analysis.Guyonnet15), 17
fwcentroid() (in module vison.point.lib), 32

G

Gaussmeter (class in vison.point.gauss), 29
gen_render() (vison.eyegore.eyegore.ImageDisplay method), 26
generate_Explog_FOCUS00() (in module vison.point.FOCUS00), 41
generate_Fake_FOCUS00() (in module vison.point.FOCUS00), 42
generate_FITS_FOCUS00() (in module vison.point.FOCUS00), 42
generate_GaussPSF() (in module vison.analysis.Guyonnet15), 17
generate_header() (in module vison.support.latex), 49
generate_HK_FOCUS00() (in module vison.point.FOCUS00), 42
generate_Latex() (vison.support.report.Figure method), 48, 50
generate_Latex() (vison.support.report.Section method), 48, 50
generate_Latex() (vison.support.report.Table method), 49, 50
generate_Latex() (vison.support.report.Text method), 49, 50
genExpLog() (in module vison.scripts.vis_genDataSet), 45
get_basic_spot_FOCUS00() (in module vison.point.FOCUS00), 42
get_centroid() (vison.point.photom.Photometer method), 30
get_cutout() (vison.datamodel.ccd.CCD method), 11

get_data() (vison.eyegore.eyegore.ExpLogDisplay method), 26
get_data() (vison.eyegore.eyegore.HKDisplay method), 26
get_data() (vison.eyegore.eyegore.ImageDisplay method), 26
get_deltaQ() (in module vison.analysis.Guyonnet15), 17
get_FOCUS00_structure() (in module vison.point.FOCUS00), 42
get_FW_ID() (in module vison.ogse.ogse), 27
get_Gauss_spot() (in module vison.point.lib), 33
get_ilum() (in module vison.pipe.FlatFielding), 10
get_ilum_splines() (in module vison.pipe.FlatFielding), 10
get_kernel() (in module vison.analysis.Guyonnet15), 17
get_mask() (vison.datamodel.ccd.CCD method), 11
get_photom_spot() (in module vison.point.lib), 33
get_quad() (vison.datamodel.ccd.CCD method), 11
get_Quadrupoles_spot() (in module vison.point.lib), 33
get_Rdisp() (in module vison.analysis.Guyonnet15), 17
get_shape_spot_FOCUS00() (in module vison.point.FOCUS00), 42
get_stats() (vison.datamodel.ccd.CCD method), 12
getacrosscolscut() (in module vison.datamodel.QLAtools), 14
getacrossrowscut() (in module vison.datamodel.QLAtools), 14
getsectioncollims() (vison.datamodel.ccd.CCD method), 12
getsectionstats() (in module vison.datamodel.QLAtools), 14
grab_numbers_and_codes() (in module vison.support.ET), 47

H

HKDisplay (class in vison.eyegore.eyegore), 26
HKplot() (in module vison.datamodel.HKtools), 13

I

ImageDisplay (class in vison.eyegore.eyegore), 26
iniExplog() (in module vison.datamodel.EXPLOGtools), 12
iniExplog() (vison.datamodel.EXPLOGtools.ExpLogClass method), 12
iniHK_QFM() (in module vison.datamodel.HKtools), 13
isNumeric() (vison.eyegore.eyegore.ExpLogDisplay method), 26
isOpen() (vison.support.ds9.ds9class method), 47

L

launch() (vison.support.ds9.ds9class method), 47
launchtask() (vison.pipe.master.Pipe method), 9
loadExpLog() (in module vison.datamodel.EXPLOGtools), 12

loadHK_preQM() (in module vison.datamodel.HKtools),
13

loadHK_QFM() (in module vison.datamodel.HKtools),
13

M

measure_bgd() (vison.point.photom.Photometer method),
30

measureRefinedEllipticity() (vison.point.shape.Shapemeter method), 31

mergeExpLogs() (in module vison.datamodel.EXPLOGtools), 12

meta_analysis() (in module vison.dark.BIAS01), 36

meta_analysis() (in module vison.dark.DARK01), 37

meta_analysis_FOCUS00() (in module vison.point.FOCUS00), 42

meta_analysis_PSF0X() (in module vison.flat.PTC0X),
40

meta_analysis_PSF0X() (in module vison.point.PSF0X),
43

P

parse_fits() (vison.pipe.FlatFielding.FlatField method),
10

parseDTstr() (in module vison.datamodel.HKtools), 13

parseHKfiles() (in module vison.datamodel.HKtools), 13

parseHKfname() (in module vison.datamodel.HKtools),
13

Photometer (class in vison.point.photom), 30

Pipe (class in vison.pipe.master), 9

plot_map() (in module vison.analysis.Guyonnet15), 17

plot_maps_ftheta() (in module vison.analysis.Guyonnet15), 17

plotAcCOLcuts() (in module vison.datamodel.QLAtools), 14

plotAcROWcuts() (in module vison.datamodel.QLAtools), 14

plotQuads() (in module vison.datamodel.QLAtools), 14

prep_data() (in module vison.dark.BIAS01), 36

prep_data() (in module vison.dark.DARK01), 37

prep_data_FOCUS00() (in module vison.point.FOCUS00), 42

prep_data_NL01() (in module vison.flat.NL01), 39

prep_data_PSF0X() (in module vison.point.PSF0X), 43

prep_data_PTC0X() (in module vison.flat.PTC0X), 40

produce_IndivFlats() (in module vison.pipe.FlatFielding),
10

produce_MasterFlat() (in module vison.pipe.FlatFielding), 10

produce_SingleFlatfield() (in module vison.pipe.FlatFielding), 10

Q

quadrupoles() (vison.point.shape.Shapemeter method),

32

R

reportFITS() (in module vison.datamodel.QLAtools), 14

reportHK() (in module vison.datamodel.HKtools), 14

run() (in module vison.flat.NL01), 39

run() (in module vison.flat.PTC0X), 40

run() (in module vison.point.FOCUS00), 42

run() (in module vison.point.PSF0X), 43

run() (vison.pipe.master.Pipe method), 9

S

search_EXPLOG() (vison.eyegore.eyegore.ExpLogDisplay method),
26

search_HKfile() (vison.eyegore.eyegore.HKDisplay method), 26

Section (class in vison.support.report), 48, 50

select_HKkeys() (vison.eyegore.eyegore.HKDisplay method), 26

set_quad() (vison.datamodel.ccd.CCD method), 12

setUpLogger() (in module vison.support.logger), 49

Shapemeter (class in vison.point.shape), 31

show_disps_CCD273() (in module vison.analysis.Guyonnet15), 17

show_spots_allCCDs() (in module vison.point.display),
29

sim_window() (vison.datamodel.ccd.CCD method), 12

simadd_flatilum() (vison.datamodel.ccd.CCD method),
12

simadd_points() (vison.datamodel.ccd.CCD method), 12

simadd_poisson() (vison.datamodel.ccd.CCD method),
12

simadd_ron() (vison.datamodel.ccd.CCD method), 12

SimpleLogger (class in vison.support.logger), 49

solve_for_A_linalg() (in module vison.analysis.Guyonnet15), 17

solve_for_psmooth() (in module vison.analysis.Guyonnet15), 18

sortBy() (vison.eyegore.eyegore.ExpLogDisplay method), 26

Spot (class in vison.point.spot), 32

SpotBase (class in vison.point.basis), 29

sub_bgd() (vison.point.photom.Photometer method), 30

sub_bias() (vison.datamodel.ccd.CCD method), 12

sub_offset() (vison.datamodel.ccd.CCD method), 12

synthHK() (in module vison.datamodel.HKtools), 14

T

Table (class in vison.support.report), 48, 50

test() (in module vison.datamodel.EXPLOGtools), 12

test() (in module vison.support.monitor), 49

test0() (in module vison.analysis.Guyonnet15), 18

test_create_from_scratch() (in module vison.datamodel.ccd), 12
 test_getkernel() (in module vison.analysis.Guyonnet15), 18
 test_load_ELVIS_fits() (in module vison.datamodel.ccd), 12
 test_selfconsist() (in module vison.analysis.Guyonnet15), 18
 test_solve() (in module vison.analysis.Guyonnet15), 18
 Text (class in vison.support.report), 49, 50

V

vison.analysis.ellipse (module), 15
 vison.analysis.Guyonnet15 (module), 15
 vison.dark.BIAS01 (module), 36
 vison.dark.DARK01 (module), 37
 vison.datamodel.ccd (module), 11
 vison.datamodel.EXPLOGtools (module), 12
 vison.datamodel.HKtools (module), 13
 vison.datamodel.QLAtools (module), 14
 vison.eyegore.eyegore (module), 25
 vison.flat.FLAT0X (module), 37
 vison.flat.NL01 (module), 38
 vison.flat.ptc (module), 21
 vison.flat.PTC0X (module), 39
 vison.image.bits (module), 23
 vison.inject.CHINJ01 (module), 35
 vison.inject.CHINJ02 (module), 35
 vison.inject.lib (module), 19
 vison.inject.plot (module), 19
 vison.ogse.ogse (module), 27
 vison.other.PERSIST01 (module), 40
 vison.pipe.FlatFielding (module), 10
 vison.pipe.master (module), 9
 vison.point.basis (module), 29
 vison.point.display (module), 29
 vison.point.FOCUS00 (module), 40
 vison.point.gauss (module), 29
 vison.point.lib (module), 32
 vison.point.models (module), 30
 vison.point.photom (module), 30
 vison.point.PSF0X (module), 42
 vison.point.shape (module), 31
 vison.point.spot (module), 32
 vison.pump.TP01 (module), 44
 vison.pump.TP02 (module), 44
 vison.scripts.HKmonitor (module), 45
 vison.scripts.quickds9 (module), 45
 vison.scripts.vis_genDataSet (module), 45
 vison.scripts.vis_mkscripts (module), 46
 vison.support.ds9 (module), 47
 vison.support.ET (module), 47
 vison.support.files (module), 48
 vison.support.latex (module), 49

vison.support.logger (module), 49
 vison.support.monitor (module), 49
 vison.support.report (module), 48, 49

W

wait_and_run() (vison.pipe.master.Pipe method), 9
 write() (vison.support.logger.SimpleLogger method), 49
 writeFITS() (vison.point.shape.Shapemeter method), 32
 writeto() (vison.datamodel.ccd.CCD method), 12
 writeto() (vison.datamodel.EXPLOGtools.ExpLogClass method), 12

X

xpaset() (vison.support.ds9.ds9class method), 47
 xpaset() (vison.support.ds9.ds9class method), 47

Z

zoomhere() (vison.support.ds9.ds9class method), 47